

PORTO-DENT II

Mobile Evacuation Unit

Compact

Measures only 13" x 16" x 26". It is no obstacle at chair-side and stores inconspicuously. Shipping weight – 40 lbs.

Value Priced

Easily the least expensive of oral evacuation units. Excellent for limited budgets.

Ideal for Backup

A lifesaver when the central system breaks down. The user just rolls it out, plugs it in, and keeps operating.

Easy Maintenance

Non-corrosive collection canister cleans easily and is equipped with an automatic overflow shut-off. Motor comes with a one-year guarantee and can be replaced in minutes with a screwdriver and wrench.

Efficient

Adjustable motor speed can evacuate indefinitely at the rate of 75 cfm, or can be adjusted to the users preference.

Special Uses

Needed especially by the dentist in a temporary location or when it is not feasible to install a central system. Perfect for such specialists as periodontists, orthodontists, endodontists, hygienists and other professions.

Moves Easily

Rolls quietly on large polyethylene carpet casters.

Cabinet

White Granite impact resistant polyethylene plastic

Specifications

By-pass motor 7.5 Amps – 115V 60 Hz. Also available in 220/230V.



Includes two hoses and two mouthpiece adaptors.

Porto-Dent is perfect where central installation is impractical, and for practitioners with an extra office to furnish.

STERI-DENT


www.sterident.com

Porto-Dent II




101008 (PV-5)-115V
101009 230V

Central VAC w/
Brushless Motor

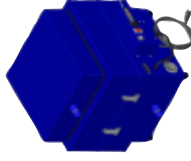


400616-115V
400617-230V

Central VAC

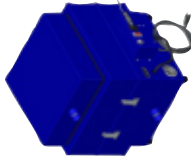


400016 (SV-2)-115V
400016-230V

Central VAC w/
Canister & Brushless Motor

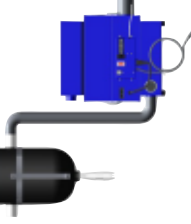

400613-115V
400614-230V

Central VAC w/ Canister

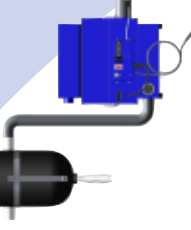

400018 (SV-4) 115V
400018-230V

SPECIFICATIONS

Electrical	7/8 HP, 8 Amps 115 VAC, 50/60 Hz 230 VAC, 50/60 Hz 7.5 Amps	1-1/2HP, 14 Amps 2-stage brushless motor 115VAC, 50/60 Hz 230 VAC, 50/60 Hz 7.2 Amps	7/8 HP, 8 Amps 115 VAC, 50/60 Hz 230 VAC, 50/60 Hz 7.5 Amps	1-1/2 HP, 14 Amps 2-stage brushless motor 115 VAC, 50/60 Hz 230 VAC, 50/60 Hz 7.2 Amps	7/8 HP, 8 Amps 115 VAC, 50/60 Hz 230 VAC, 50/60 Hz 7.5 Amps	400018 (SV-4) 115V 400018-230V
Motor Control	Variable Speed Control	Remote Motor Start	Remote Motor Start	Remote Motor Start	Remote Motor Start	Remote Motor Start
Construction	Molded polyethylene Impact Resistant	Sheet Metal Powder Coated	Sheet Metal Powder Coated	Sheet Metal Powder Coated	Sheet Metal Powder Coated	Sheet Metal Powder Coated
Chair Capacity	Comes with 2 hoses	2-5	1-2	2-5	1-2	1-2
Vacuum Suction	80" of water (6" of Hg)	120" of water (9" of Hg)	80" of water (6" of Hg)	120" of water (9" of Hg)	80" of water (6" of Hg)	80" of water (6" of Hg)
Vacuum Airflow	75 cfm	125 cfm	75 cfm	125 cfm	75 cfm	75 cfm
Canister	Stainless Steel -2 Liter	None	None	Molded plastic - 5L	Molded plastic - 5L	Molded plastic - 5L
Weight (unit+canister)	35 lbs	52 lbs	46 lbs	60 (52 + 8) lbs	54 (46 + 8) lbs	54 (46 + 8) lbs
Shipping Weight	40 lbs	70 lbs	70 lbs 2 pkgs	100 lbs 2 pkgs	100 lbs 2 pkgs	100 lbs 2 pkgs
Ships via	UPS	UPS	UPS	UPS	UPS	UPS
Dimensions (unit)	13"W x 16"D x 26"H	17"W x 19"D x 19"H	17"W x 19"D x 19"H	17"W x 19"D x 19"H	17"W x 19"D x 19"H	17"W x 19"D x 19"H
Dimensions (canister)	7" Dia. x 8"H	N/A	N/A	10" Dia. x 16"H	10" Dia. x 16"H	10" Dia. x 16"H
Shipping Dim.	16" W x 20" D x 30" H	25" W x 25" D x 25" H	25" W x 25" D x 25" H 8" W x 8" O x 8" H	25" W x 25" D x 25" H 30" W x 15" D x 15" H	25" W x 25" D x 25" H 30" W x 15" D x 15" H	25" W x 25" D x 25" H 30" W x 15" D x 15" H
Options	Portable amalgam separator Carry Handles Heavy-Duty Casters External Vent Kit for Nitrous Oxide	Amalgam separator	Amalgam separator	Amalgam separator	Amalgam separator	Amalgam separator

Note: All units available in 220-240 VAC, 50/60Hz