


Cox Dry Heat Sterilizers

Sterilization Cycles in as little as 6 Minutes

- **Fast Turn Around**
- **Eliminate Corrosion or Rusting of Instruments**
- **Low Maintenance**
- **No Supplies Needed**

Benefits of Dry Heat Sterilization

- The COX RAPID DRY HEAT STERILIZER utilizes dry heat rather than steam and pressure to sterilize. It won't corrode, dull or pit your instruments or cause them to rust, thus extending the life of every costly instrument.
- No need to purchase supplies such as distilled water, rust inhibitors, or chemical solutions.
- The superior speed of the COX STERILIZER gets instruments turned around to you quicker, which means avoiding the cost of buying extra sets of instruments.

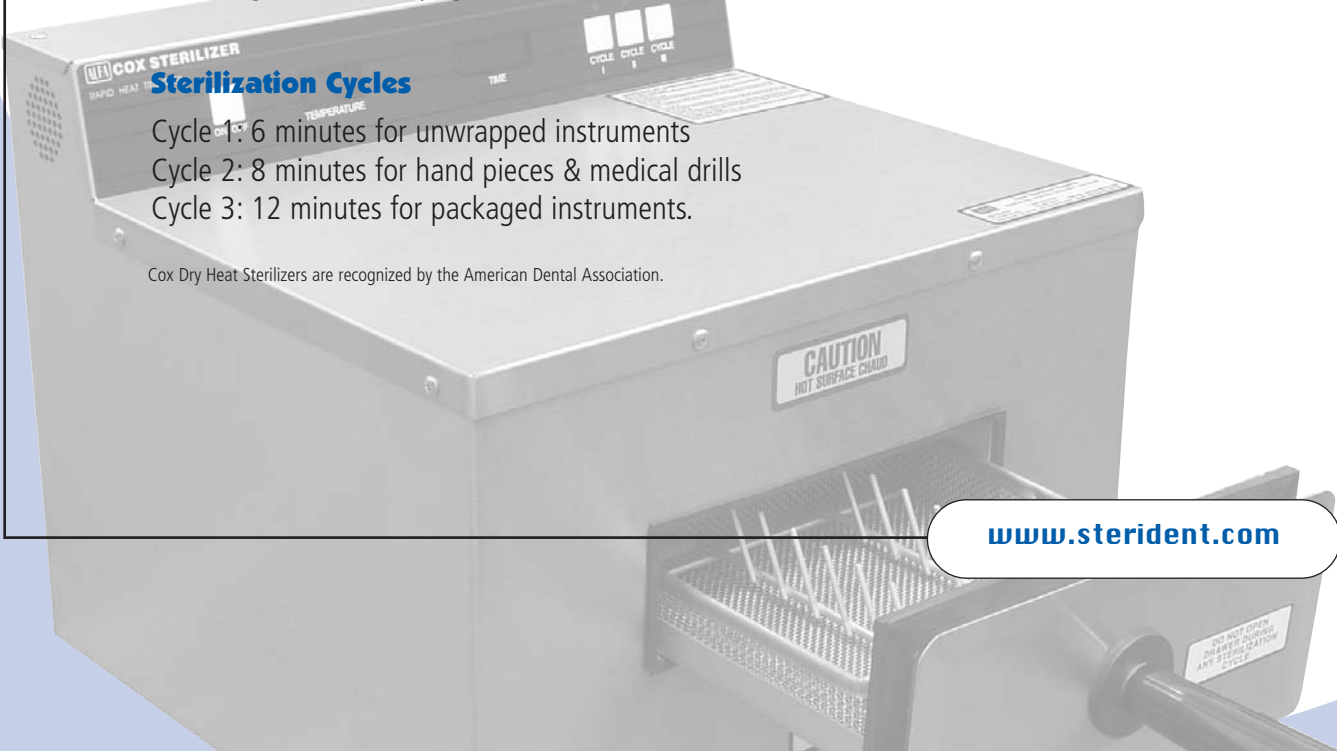

Sterilization Cycles

- Cycle 1: 6 minutes for unwrapped instruments
- Cycle 2: 8 minutes for hand pieces & medical drills
- Cycle 3: 12 minutes for packaged instruments.

Cox Dry Heat Sterilizers are recognized by the American Dental Association.

www.sterident.com

Fast, Easy to Use and Convenient

Specifications

Operating Cycles:	6 minute, 8 minute, or 12 minute
Unit Size:	17"W x 17"D x 15"H 432mm x 432mm x 381mm high
Weight:	50 lbs. (23Kg)
Instrument Basket Dimensions:	7.5"D x 8.25"W x 1.5"H 229mm x 216mm x 38mm
Electrical:	110V, 60Hz, 1100 Watts warm up, 300 Watts operating 220V, 50Hz, 1100 Watts warm up, 300 Watts operating
Operating Temperature:	375° Fahrenheit (190° Celsius)
Warranty:	One year (Parts)
Models:	COX-115V-N COX-220V-N
Approvals:	UL Listed 510 (k) K881371

Standard Accessories:

One Drawer
One Removable Instrument Basket
Basket Removal Tool
Cooling Rack
Instrument Rack (for packaged instruments)
Mesh Burr Holder

Optional:
Handpiece Processing Kit
(P/N CS6030-A)
External Cooling System
(P/N PCOX0001)

cpac
Imaging

Toll Free: (800) 828-6011
Tel: (585) 382-3223
Fax: (585) 382-9481